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
## Storage information

The storage location should be dry and at room temperature. The product has to be protected against direct sunlight exposure and humidity.

## Cleaning information

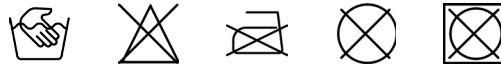Wash by hand using gentle detergents in temperature max 35 ° C. Do not chlorinate, do not dry clean, wring gently before drying. Do not dry close to the artificial warmth source. Do not iron.

## Guarantee

The guarantee lies within the framework of the legal regulations.

## Disposal

The product can be disposed of with the normal household waste.

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 mdh sp. z o.o.  
ul. Maratońska 104  
94-007 Łódź  
www.qmedbymeira.pl  
ortopedia@mdh.pl  
tel. +48 42 636 81 44

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(EN)

### Acromed shoulder brace

Thank you very much for selecting a product of the Meyra-Ortopedia Group. The product complies with all legal regulations according to the Medical Products Law.

Please read the directions for use carefully and keep to the advice and recommendations given. Should you have further questions please address the responsible physician, your local specialist dealer or our service hotline.

### Product description

Breathable fabric provides optimal comfort of use. Proper compression of soft tissues protects against oedema and reduces muscle tension. System of bands ensures stabilization of shoulder and acromioclavicular joints, enables adjustment to individual needs. Additional, flexible bands allow to limit internal and external rotation. Product is available in two versions: left and right.

### Safety instructions

The first application of the product must be carried out by special personnel trained on these kind of products as otherwise the product liability cannot be guaranteed. The product must not be worn on open wounds.

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## Indications:

Post-operative treatment of humeral head fracture, glenoid cavity, shoulder joint dislocation  
Degeneration of shoulder joint and related pain syndromes.  
Habitual dislocation of shoulder joint  
Follow-up treatment after reconstruction of labrum glenoidale and articular capsule.  
Protection against shoulder joint dislocation in case of upper limb paralysis of various etiology

## Directions for use

Always remember that choosing proper size of the support is essential for its effectiveness.  
Pull the support so that the head of humerus would be placed in designed for it area.

Fasten the inflexible bands around the chest and adjust the support to the chest circumference.  
Then, fasten the flexible bands according to therapeutic need:

- If you fasten the front and rear bands symmetrically it will limit abduction range of motion.
- If you fasten the front band stronger, it will limit external rotation.
- If you fasten the rear band stronger, it will limit internal rotation.

## Contraindications

Conditions which require shoulder joint immobilization  
Wounds and skin changes

## Wear time

Wear time is determined in consultation with your doctor.

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## Mechanism of acting

Construction ensures correct positioning of humeral bone. Properly placed compression area protect against oedema, provide relief of muscles and ligaments. Additional bands limit abduction and internal or external rotation (depending on the strength and way of fastening).

## ACROMED

shoulder support effectively helps in treatment and rehabilitation. It allows to position properly the head of humerus in the acetabulum, to relieve the muscles, ligaments of shoulder and acromioclavicular joints. It also enables limitation of abduction range of motion.

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